AHP2

Indicator Strips

For use with 2% High Level
Disinfectant and Sterilant

Part No. AHP2

Description

Accelerated Hydrogen Peroxide® (AHP®) Indicator Test Strips are easy to use dip-and-read reagent strips for a pass or fail determination of hydrogen peroxide concentration in disinfectant solutions for variety of applications.

Direction For Use

1. Dip the padded end of the indicator test strips into the test solution and remove after one second.
2. Gently shake off the excess solution. Once removed, hold the strip horizontally and start the timer.
3. Wait for 60 seconds then compare the colour on the reagent pad with the reference colour block. If the colour on the reagent pad is lighter than the reference colour block for a "PASS", then the concentration of hydrogen peroxide is below the minimum recommended concentration for the product and should be discarded. "FAIL" colour blocks are provided only for reference and should not be used for measuring or quantifying the level of hydrogen peroxide in the solution.

NOTE: Delayed reading of the strip color may invalidate the results.

Shelf life

Do not use the strip after the expiration date. Please refer to expiry date on indicator test strip bottle.

Storage

Store the Accelerated Hydrogen Peroxide® Indicator Test Strips at room temperature between 59°F and 86°F (15°C–30°C). Do not expose the strip to direct sun light. Ensure the container is closed when not used.

Important Precautions

Keep the strip in the original container. Remove only the required strip before testing and immediately replace the cap tightly. Exposure of the uncovered strip to light and moisture may cause discoloration of the strip. Do not touch the reagent pad with finger or any other surface.

Certificate of Compliance is available upon request.

Principle of the Test

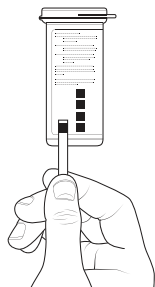
The reaction of the AHP® Test is based on a proprietary coupling system where the product reacts with an iodide salt and a coupler to form a brown coloured complex.

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Accelerated Hydrogen Peroxide® and Design are trademarks of Diversey, Inc.
Virox Technologies Inc.