Indicator Strips

Part No. AHP425

Indicator Strips For Use with 4.25%
Concentrate AHP® Products
(1:16, 1:32, 1:64, 1:128 Dilution)

Description

Accelerated Hydrogen Peroxide® (AHP®) Indicator Test Strips are easy to use, dip-and-read reagent strips to verify the concentration of hydrogen peroxide, the active ingredient in AHP®. **The Test Strips were DESIGNED for use with the 4.25% AHP® Surface Disinfectant Concentrates diluted to 1:16 1:32, 1:64, and 1:128.**

Direction For Use

(1:16, 1:32, 1:64, 1:128)

1. Dip the padded end of the indicator test strips into the diluted solution and remove after one second.
2. Do not shake the test strip. Once removed, hold the strip horizontally and start the timer.
3. Match the indicator test strip to one of the following color blocks between 50 and 55 seconds to verify the dilution rate.

NOTE: Delayed reading of the strip color may invalidate the results.
The strip color will continue to darken after the read time.

Shelf life

Do not use the strip after the expiration date (opened or unopened).
Please refer to expiry date on indicator test strip bottle.

Storage

Store the Accelerated Hydrogen Peroxide® (AHP®) indicator test strips at room temperature between 15°C–30°C (59°F–86°F). Do not expose the strip to direct sunlight. Ensure the cap is closed tightly after each use.

Important Precautions

Keep the unused test strips in the original container. Remove only the required test strip before testing and ensure the cap is closed tightly after each use. Exposure of the uncovered test strip to light and moisture may cause discoloration of the strip. Do not touch the reagent pad with finger or any other surface.

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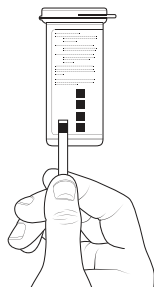
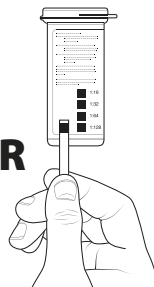
REF274 (26/264)

Accelerated Hydrogen Peroxide® and Design are trademarks of Diversey, Inc.

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Virox Technologies Inc.